



Gill III

MODULAR ADJUSTABLE SEATING SYSTEM

Standard Gill III comprises of:

- Headrest
- Pelvic Belt
- Seat board c/w Lock & Latch interface
- 15° ramped cushion
- Backrest c/w Lock & Latch interface
- Fixed Pelvic Supports
- Swing away thoracic supports

SIZE GUIDE	STANDARD SIZES LISTED BELOW								
Gill III model	A	B	C	D	E	F	G	H	I
Wheelchair size	305 (12")	330 (13")	355 (14")	380 (15")	405 (16")	430 (17")	455 (18")	480 (19")	510 (20")
Seat width min	180	190	215	240	265	290	315	340	365
Seat width max	280	290	315	340	365	390	415	440	465
Seat depth min	220	245	245	325	325	350	375	375	375
Seat depth max	345	360	360	440	440	465	490	490	490
Backrest height	360	410	410	460	460	510	510	560	560

Please note: the above dimensions are in conjunction with 8 & 9 series wheelchairs. These may vary due to other wheelchair manufacturers tolerances. All dimensions are in mm unless otherwise stated

Available in seven colours



- **Anti-bacterial/anti-fungal:** will not support microbial or fungal growth including Salmonella, E Coli and MRSA (Methicillin Resistant Staphylococcus Aureus)
- **Waterproof:** Water or oil based solutions will not penetrate or absorb into the fabric. Hydrostatic Head BS3424 > 1 metre
- **Breathable:** Moisture vapour permeable BS3546 65%
- **Easy to maintain:** Extremely Stain Resistant. Spillages can be easily removed using a damp sponge or cloth and wiping dry
- **Flame retardant:** BS5852 Crib 5
- **Stain resistant:** Excellent stain removal has been observed for blood, urine, grease, ink, coffee, iodine, betadine, ketchup, chewing gum, chocolate, grape juice

Using the latest materials, the Gill III seating system can be fully reconditioned.

See our Accessory Catalogue for the full range of belts, harnesses, headrests, etc. available.



Successfully crash tested to ISO 16840 part 4

Please see our current price list of the full range of standard equipment and accessories. All rights reserved to design changes.



Gill III SYSTEM



EDITION 2



Rehabilitation Manufacturing Services Limited



Medway House, 277 Gillingham Road, Gillingham, Kent, ME7 4QX
 Tel: +44 (0)1634 578881 Fax: +44 (0)1634 581456
 Email: sales@rms-kent.co.uk www.ineedawheelchair.co.uk
 Rehabilitation Engineering For a Caring World



Gill III


MODULAR ADJUSTABLE SEATING SYSTEM

A completely versatile way to suit your clients needs whether you prescribe the complete system or individual components.

Successfully crash tested to ISO 16840 part 4

Gill III - Seat Base

- Provides a lightweight and firm seat base
- Manufactured from black 12mm polypropylene
- Interfaced into wheelchair using 15mm drop in brackets and located using Lock & Latch interface system
- Fixing holes located along entire length to enable seat depth adjustment and interfacing
- Location holes for pommel or kneeblock mounting

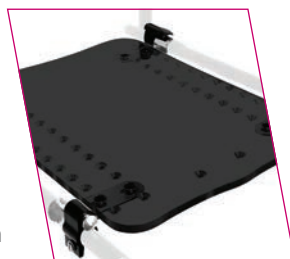

Gill III - Standard Cushion

- Available flat or with 10°/ 15° ramp
- Platilon inter-liner with removable covers
- RG60 foam
- Pressure relief foams available



Gill III - Kneeblocks/ Pommel Attachment

- Attached to seat base by either a ratchet or pop-pin mounting
- Available in various sizes
- Adjustable in width, height and depth
- Can accommodate asymmetrical leg lengths



Kneeblock



Pommel


Gill III - Headrest

- Moulded headrest which gives basic support
- Height, depth and angle adjustments as standard
- Alternative headrests available in RMS Accessories catalogue



Gill III - Backrest

- Flat backrest manufactured from black 12mm polypropylene
- Backrest interfaced using Lock & Latch system
- Interface brackets have 75mm depth and 25° backrest angle adjustment
- Backrest cushion has Platilon inter-liners with removable covers

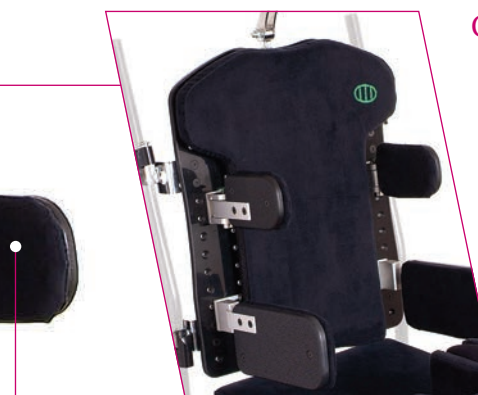

Gill III - Thoracic Pads

- Available in either swing away or fixed
- Choice of pad sizes
- Available flat or curved
- Easily adjusted
- 100mm of width adjustment



Gill III - Pelvic Supports

- Available in either swing away or fixed
- Choice of pad sizes
- Easily adjusted
- 100mm of width adjustment



For more information call: 01634 578881