

DORY

TILT-IN-SPACE
WHEELCHAIR



DORY TILT-IN-SPACE WHEELCHAIR

The Dory is a tilt-in-space wheelchair base with a lightweight aluminium frame that is the perfect wheelchair base to suit any situation. Compatible with virtually any specialist seating system; with seat extension rails to accommodate various seat depths; and with a massive range of 50° tilt-in-space, this makes the Dory one of the most versatile wheelchairs on the market.

KEY FEATURES

- Ideal wheelchair for variety of applications
- Choose between a folding or fixed frame
- Three sizes available with width adjustment
- Up to 55cm seat depth and 50° tilt-in-space
- Adjustable centre of gravity for stability
- Fully transport crash tested to ISO-7176-19



ABOUT THE DORY

The Dory comes with a great selection of equipment as standard. A black folding or solid aluminium frame that is width and depth adjustable for growth; height adjustable handles; 300mm solid rear wheels and 150mm front castor wheels; push lock brakes; and tie down points for transportation. Alongside this there are many upgrades available including various wheel sizes; armrests and footrests; RMS manufactured belts and harnesses; and trays and other accessories.





TILTING FRAME

The Dory offers a extensive range of movement in the frame with 50° of tilt-in-space. The frame can be ordered in two standard configurations, one with 0° to 50° range and other with -5° to 45° or range allowing for forward tipping - making it easier for users to transfer from the wheelchair.

FOLDING OR FIXED

The Dory comes available with either a rigid or folding frame - this means the Dory is perfect for a wide range of clinical settings - either a more secure option to be permanently set up, or the folding option which can be easily folded and packed away for transportation.




VERSATILE INTERFACE

The Dory frame is extremely adaptable wheelchair frame that can be fixed to virtually any third party seating system. The Dory is even better suited to match with our **Caddy Modular Seating System** or the **Virtus Dynamic Seating System** - contact us to find out more about these products.

STANDARD PARTS

- Black aluminium frame
- 50° tilt-in-space
- Height adjustable push handles
- 150mm solid front castors
- 300mm solid rear wheels
- Push lock brakes
- Tie down points for transportation

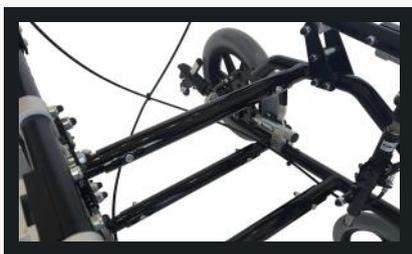

The Dory has a width adjustable base (**shown here on the left**) that allows the wheelchair base to be widened or narrowed for great flexibility.

The adjustment is simple; easy to make; and allows fantastic growth in the frame.

EXTRA PARTS

- Various wheel + castor sizes
- Hub brakes
- Armrests
- Footrests
- Moulded tray
- Range of belts + harnesses
- Ventilator + oxygen carriers

TECHNICAL SPECIFICATIONS

	Small	Medium	Large
Seat width	31 - 37	37 - 43	43 - 50
Seat depth	35 - 45	40 - 50	45 - 55
Seat angle	0° to 50° or -5° to 45°		
Seat rail height	43 - 45 - 46 - 48		
Backrest angle	90° - 130° (4 positions)		
Maximum user weight	135 kg		
Total chair weight	18 - 25 kg		



Thompson House, Unit 10, Styles Close,
Sittingbourne, Kent, ME10 3BF.

Tel: +44 (0) 1795 477280
Fax: +44 (0) 1795 229692
Email: sales@rms-kent.co.uk

www.rms-rehab.co.uk

   @RMSRehab